

OSSTEM[®]
IMPLANT

2020 - 21
OS SYSTEM
CATALOG

CE PG
2460 ИМ22 ENG_ver.

OSSTEM[®]
IMPLANT

Osstem Implant
2020-21 Comprehensive Catalog

Planning/Editing Design Center Brand Design Team

Supervision PM Headquarters, Implant Lab

Production/Distribution PM Headquarters

Date of Publication 02.2020


Publisher OSSTEM IMPLANT Co., Ltd. 3, Magokjungang 12-ro,
Gangseo-gu, Seoul 07789, (Republic of) Korea

Phone +82.2.2016.7000

Fax +82.2.2016.7001

www.osstem.com

2020 Printed in Seoul, Korea



OS SYSTEM

004 OrthAnchor Simple Head
006 OrthAnchor Through Hole
008 OrthAnchor Small Head
009 OrthAnchor Bracket Head
010 OrthAnchor Simple Head Half Etched
012 OrthAnchor Through Hole Half Etched


014 Ortho KIT
018 ORP KIT
024 e-Driver & e-Driver plus
025 V-ceph

OrthAnchor Simple Head



Simple Head

- Machined surface
- Material : Ti-6Al-4V
- Connected component : Coil spring(Ø2.5), Power chain, Elastic band
- ※ G/H 4.0 type is produced after order.

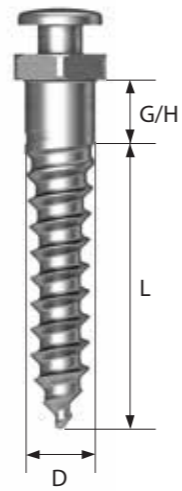

D Ø1.6

| G/H \ L | 6 | 8 | 10 |
|---------|------------|-----------|-----------|
| 1.5 | OSSH 1606 | OSSH 1608 | OSSH 1610 |
| 4.0 | OSSH 16064 | - | - |

D Ø1.2

| G/H \ L | 6 | 8 | 10 |
|---------|-----------|-----------|----|
| 1.5 | OSSH 1206 | OSSH 1208 | - |

D Ø1.8

| G/H \ L | 6 | 8 | 10 |
|---------|------------|-----------|-----------|
| 1.5 | OSSH 1806 | OSSH 1808 | OSSH 1810 |
| 4.0 | OSSH 18064 | - | - |

D Ø1.4


| G/H \ L | 6 | 8 | 10 |
|---------|-----------|-----------|----|
| 1.5 | OSSH 1406 | OSSH 1408 | - |

OrthAnchor Through Hole



Through Hole

- Machined surface
 - Material : Ti-6Al-4V
 - D (hole) : Ø0.8
 - Connected component : Arch wire(round), Coil spring(Ø2.5), Power chain, Elastic band
- ※ G/H 4.0 type is produced after order.

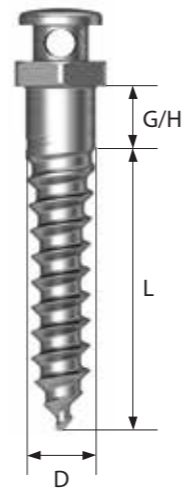

D Ø1.6

| G/H \ L | 6 | 8 | 10 |
|---------|------------|-----------|-----------|
| 1.5 | OSTH 1606 | OSTH 1608 | OSTH 1610 |
| 4.0 | OSTH 16064 | - | - |

D Ø1.2

| G/H \ L | 6 | 8 | 10 |
|---------|-----------|-----------|----|
| 1.5 | OSTH 1206 | OSTH 1208 | - |

D Ø1.4

| G/H \ L | 6 | 8 | 10 |
|---------|-----------|-----------|----|
| 1.5 | OSTH 1406 | OSTH 1408 | - |


D Ø1.8

| G/H \ L | 6 | 8 | 10 |
|---------|------------|-----------|-----------|
| 1.5 | OSTH 1806 | OSTH 1808 | OSTH 1810 |
| 4.0 | OSTH 18064 | - | - |




Small Head

- Machined surface
- Material : Ti-6Al-4V
- D (head) : $\varnothing 1.48$
- Connected component : Coil spring($\varnothing 1.5/2.0/2.5$), Power chain, Elastic band

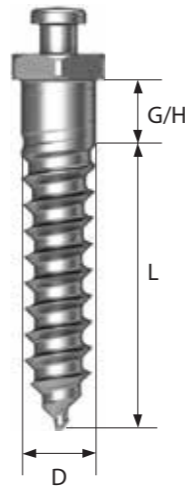

Bracket Head

- Machined surface
- Material : Ti-6Al-4V
- Excellent compatibility with various arch wires
- Easy path adjustment with the cross wire slot
- Connected component : Arch wire(rec./round), Coil spring($\varnothing 2.5$), Power chain, Elastic band

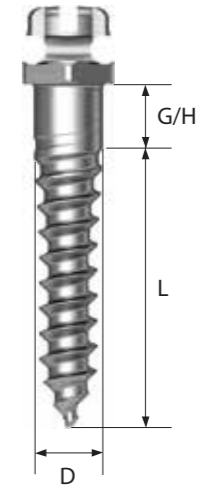

| D \varnothing | G/H \ L | 6 | 8 | 10 |
|---------------------|---------|------------|------------|------------|
| D $\varnothing 1.4$ | | | | |
| | 1.5 | OSSHs 1406 | OSSHs 1408 | - |
| D $\varnothing 1.6$ | | | | |
| | 1.5 | OSSHs 1606 | OSSHs 1608 | OSSHs 1610 |
| D $\varnothing 1.8$ | | | | |
| | 1.5 | OSSHs 1806 | OSSHs 1808 | OSSHs 1810 |


| D \varnothing | G/H \ L | 6 | 8 | 10 |
|---------------------|---------|-----------|-----------|-----------|
| D $\varnothing 1.4$ | | | | |
| | 1.5 | OSBH 1406 | OSBH 1408 | - |
| D $\varnothing 1.6$ | | | | |
| | 1.5 | OSBH 1606 | OSBH 1608 | OSBH 1610 |
| D $\varnothing 1.8$ | | | | |
| | 1.5 | OSBH 1806 | OSBH 1808 | OSBH 1810 |

OrthAnchor Simple Head Half Etched




Simple Head Half Etched


- Acid etched surface
- Material : Ti-6Al-4V
- Minimization of early drop out possibility
- Stable effect for children or adolescents or cases with poor bone quality
- Connected component : Arch wire(round), Coil spring(Ø2.5), Power chain, Elastic band

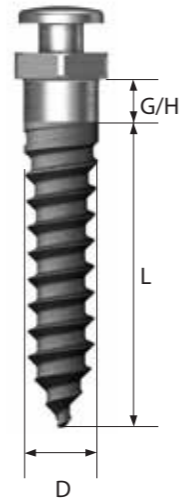

D Ø1.6

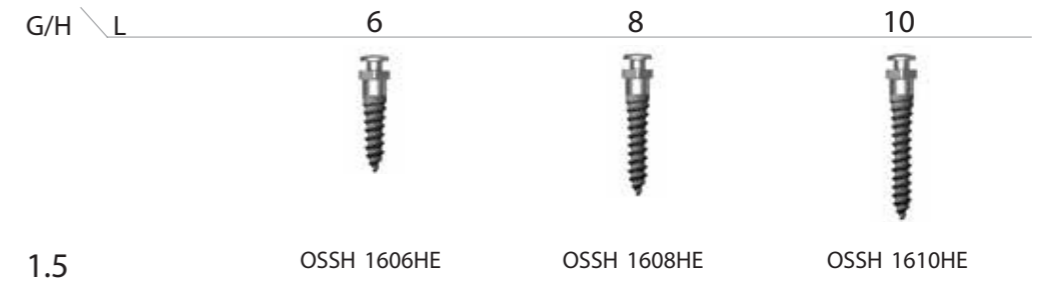

D Ø1.2

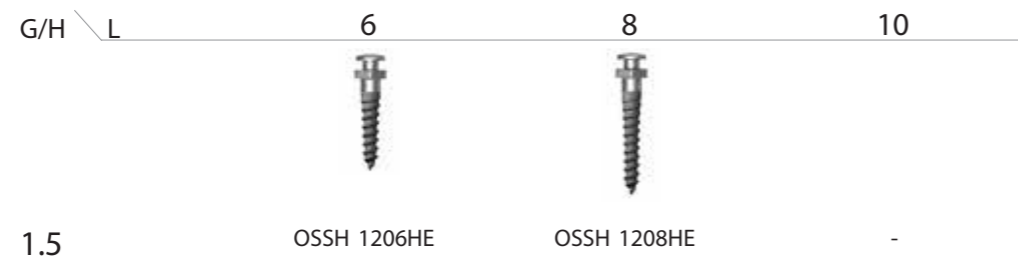

D Ø1.8



D Ø1.4




OrthAnchor Through Hole Half Etched



Through Hole Half Etched

- Acid etched surface
- Material : Ti-6Al-4V
- Minimization of early drop out possibility
- Stable effect when applying for children or adolescents or cases with poor bone quality
- Connected component : Arch wire(round), Coil spring(Ø2.5), Power chain, Elastic band

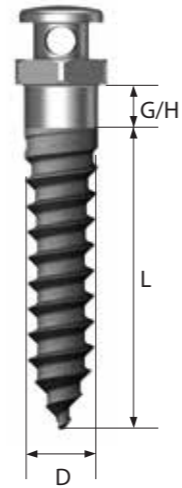

D Ø1.6

| G/H \ L | 6 | 8 | 10 |
|---------|-------------|-------------|-------------|
| 1.5 | OSTH 1606HE | OSTH 1608HE | OSTH 1610HE |

D Ø1.2

| G/H \ L | 6 | 8 | 10 |
|---------|-------------|-------------|----|
| 1.5 | OSTH 1206HE | OSTH 1208HE | - |

D Ø1.8


| G/H \ L | 6 | 8 | 10 |
|---------|-------------|-------------|-------------|
| 1.5 | OSTH 1806HE | OSTH 1808HE | OSTH 1810HE |

D Ø1.4

| G/H \ L | 6 | 8 | 10 |
|---------|-------------|-------------|----|
| 1.5 | OSTH 1406HE | OSTH 1408HE | - |




For OS

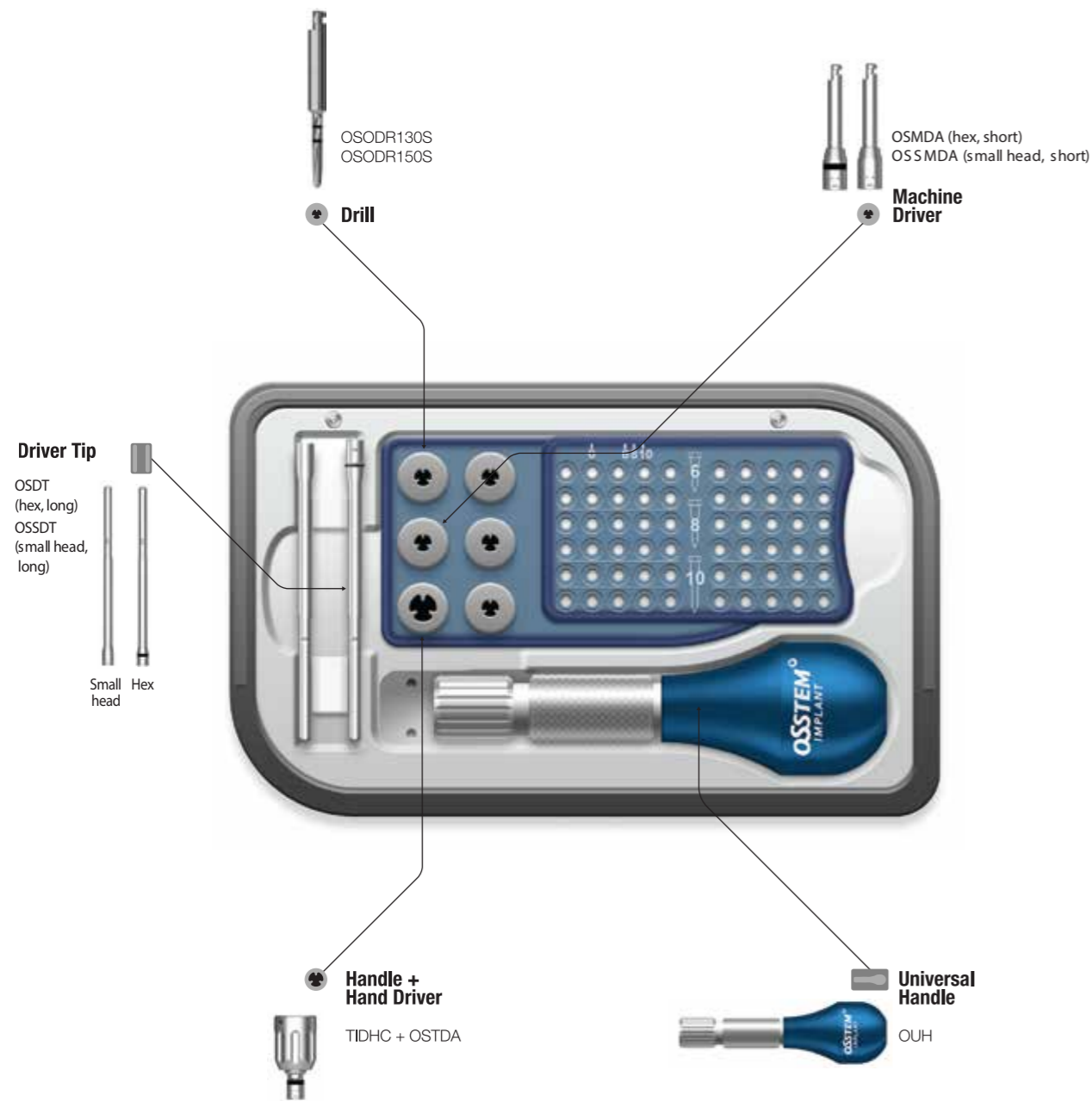

Drill 10.2013


- Connected to a hand piece (engine) for use
- Ø1.0 drill : For Ø1.2/1.4 screw procedure
- Ø1.3 drill : For Ø1.6 screw procedure
- Ø1.5 drill : For Ø1.8 screw procedure
- Recommended speed : 800rpm (high speed)
- Insertion placement recommended after removing cortical bone only
(Drilling to the same length as the screw length if the cortical bone is too thick)
- Ø1.0 drill for optional purchase (not included in the KIT)

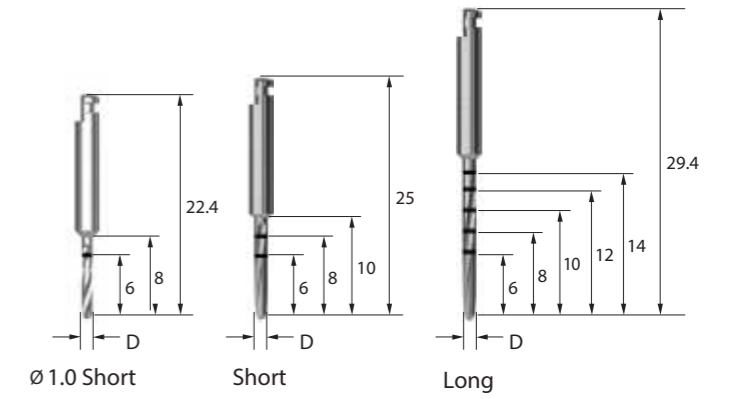

| L \ D | Ø 1.0 | Ø 1.3 | Ø 1.5 |
|-------|------------|------------|------------|
| Short | OSODR 100S | OSODR 130S | OSODR 150S |
| Long | - | OSODR 130C | OSODR 150C |

Universal Handle 01.2009


- Used after assembling to the driver tip, easy procedure with anti-slip handle

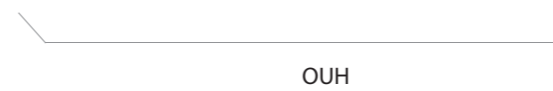

Driver Tip 01.2009

- Assembled to a universal handle for OrthAnchor procedure
- Consisted of a general hex driver and a small head driver
- Compatible with a universal handle of other company (J company)

| L \ Type | Hex | Small Head |
|-----------|------------|-------------|
| Short (L) | OSDTS (45) | OSSDTS (45) |
| Long (L) | OSDT (67) | OSSDT (67) |

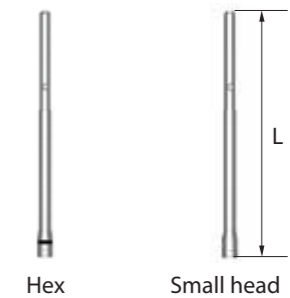

Ortho KIT Surgical Instruments

Hand Drill 08.2012

- Assembled to a universal handle for use
- Removing cortical bone only
- Drilling depth : 4mm
- Optional purchase (not included in the KIT)
- ※ Maintain the drilling direction, not exerting bending load while using



OSHDR130

Machine Driver 01.2009

- Connected to the engine for OrthAnchor surgery
- Consisted of a general hex driver and a small head hand driver

| L \ Type | Hex | Small Head |
|-----------|--------------|---------------|
| Short (L) | OSMDA (21.4) | OSSMDA (21.4) |
| Long (L) | OSMDB (31.4) | OSSMDB (31.4) |



Hex



Small head

Driver Handle 01.2009

- Used for tightening screws with a hand after connecting a hand driver



TIDHC

Hand Driver 01.2009

- Connected to a driver handle or a ratchet wrench for OrthAnchor screw procedure
- Consisted of a general hex driver and a small head hand driver
- Small head hand driver for optional purchase (not included in the KIT)




Hex




Small head

| Type | Hex | Small Head |
|------|-------|------------|
| | OSTDA | OSSTDA |

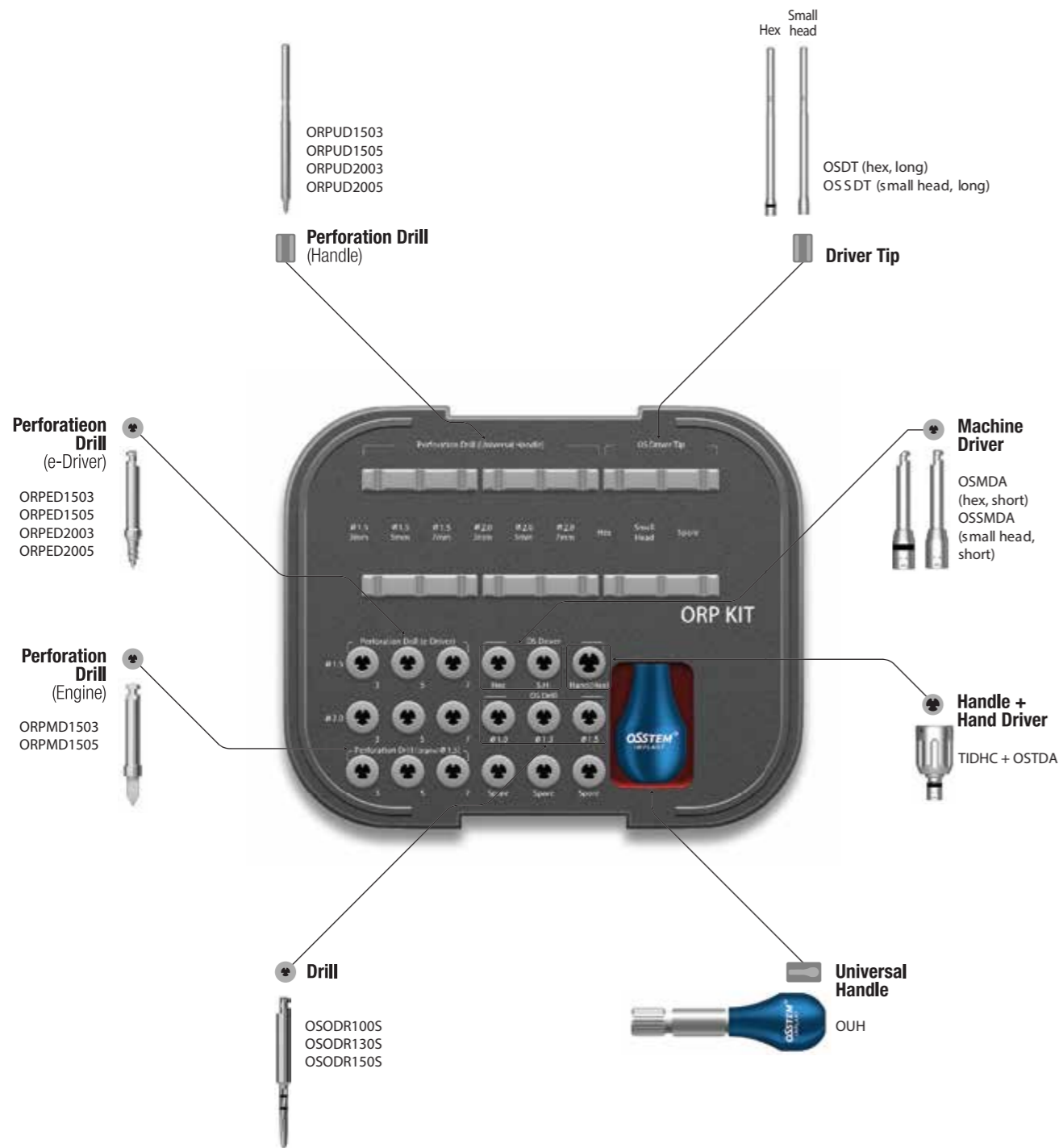

Perforation Drill (Handle)


- Connected to universal handle for MOP operation
- Used for areas easy to reach
- Ø1.5 : Anterior region / Ø2.0 : Posterior region and areas with wide spacing between teeth
- ※ MOP : micro-osteoperforation

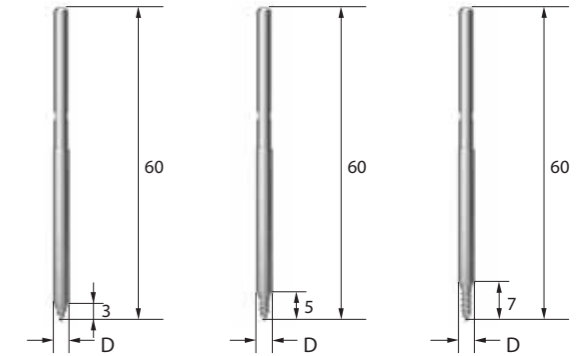

| L \ D | Ø1.5 | Ø2.0 |
|-------|------------|------------|
| 3.0 | ORPUD 1503 | ORPUD 2003 |
| 5.0 | ORPUD 1505 | ORPUD 2005 |
| 7.0 | ORPUD 1507 | ORPUD 2007 |

Perforation Drill (e-Driver)

- Connected to hand piece (engine) for MOP operation
- Used for areas hard to reach with a hand drill such as palatal region
- Recommended tightening : 25Ncm
- Recommended speed : 30-60rpm

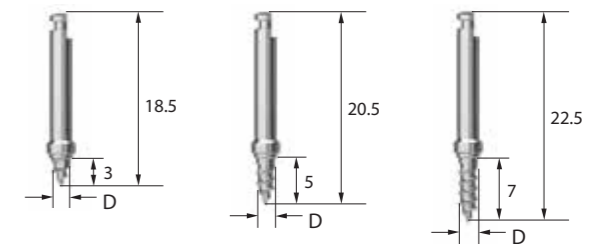


| L \ D | Ø1.5 | Ø2.0 |
|-------|------------|------------|
| 3.0 | ORPED 1503 | ORPED 2003 |
| 5.0 | ORPED 1505 | ORPED 2005 |
| 7.0 | ORPED 1507 | ORPED 2007 |

ORP KIT Surgical Instruments

Perforation Drill (Engine)

- Connected to hand piece(engine) for MOP operation (e-Driver cannot be used)
- Fast operation for hard bone or areas hard to reach with a handle drill
- Recommended speed : 1200rpm

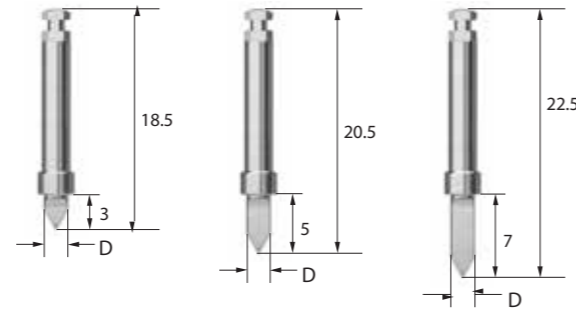

| L \ D | Ø 1.5 | Ø 2.0 |
|-------|------------|------------|
| 3.0 | ORPMD 1503 | ORPMD 2003 |
| 5.0 | ORPMD 1505 | ORPMD 2005 |
| 7.0 | ORPMD 1507 | ORPMD 2007 |

Universal Handle


- Used for MOP surgery by connecting perforation drill (for handle)
- Screw placement is available by connecting a dedicated driver tip



Driver Tip


- Used for OrthAnchor surgery by connecting universal handle
- Composed of general hex driver and small head driver
- Compatible with other company's universal handle (J, Initial of company)

| L \ Type | Hex | Small Head |
|-----------|------------|-------------|
| Short (L) | OSDTS (45) | OSSDTS (45) |
| Long (L) | OSDT (67) | OSSDT (67) |

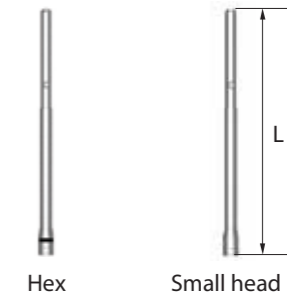

Drill

- Connected to hand piece(engine) for use
- Ø1.0 drill : used for Ø1.2/1.4 screw
- Ø1.3 drill : used for Ø1.6 screw
- Ø1.5 drill : used for Ø1.8 screw
- Recommended speed : 800rpm
- Removal of cortical bone and placement is recommended (If the cortical bone is very thick, drilling has to be same as screw length)

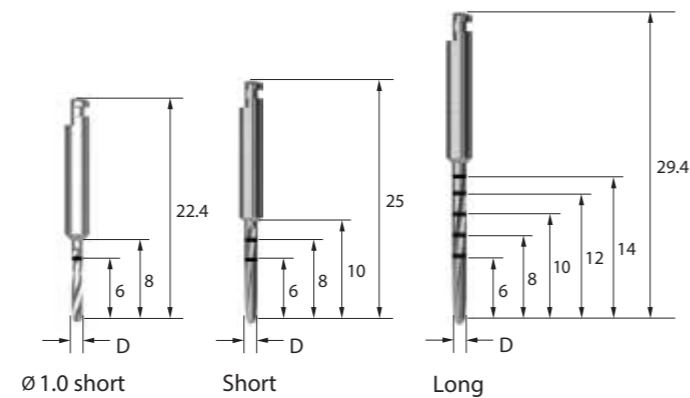

| L \ D | Ø 1.0 | Ø 1.3 | Ø 1.5 |
|-------|------------|------------|------------|
| Short | OSODR 100S | OSODR 130S | OSODR 150S |
| Long | - | OSODR 130C | OSODR 150C |

Hand Drill

- Connected to universal handle
- It can remove only cortical bone
- Drilling depth : 4mm
- Optional purchase (not included in the KIT)
- ※ Maintain drilling direction so that no bending load is applied



Driver Handle

- Use for manually fastening screws after connecting a hand driver



Hand Driver

- Used for OrthAnchor surgery by connecting to driver handle and ratchet wrench
- Composed of general hex driver and hand driver for small head
- Hand driver for small head is optional purchase (not included in the KIT)

| Type | Hex | Small Head |
|------|-------|------------|
| | OSTDA | OSSTDA |



Hex



Small head

Machine Driver

- Used in OrthAnchor operation by fastening to engine
- Composed of general hex driver and machine driver for small head

| L | Type | Hex | Small Head |
|-----------|------|--------------|---------------|
| Short (L) | | OSMDA (21.4) | OSSMDA (21.4) |
| Long (L) | | OSMDB (31.4) | OSSMDB (31.4) |




Hex



Small head

Removal Tool (Handle)

- Easily removable when fracturing OrthAnchor
- Connect to a universal handle and use it in reverse
- Select the specification according to the fracture screw diameter
- Can also be used for other company screw fractures

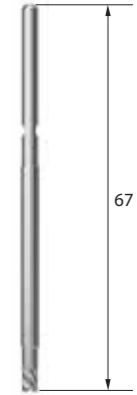

| D (Removal Screw) | Ø 1.2 | Ø 1.4 | Ø 1.6 | Ø 1.8 | Ø 2.0 |
|-------------------|----------|----------|----------|----------|----------|
| | OSRT 12H | OSRT 14H | OSRT 16H | OSRT 18H | OSRT 20H |

Removal Tool (Engine)

- Easily removable when fracturing OrthAnchor
- Connect to e-driver or hand piece (engine) and use it in reverse
- Select the specification according to the fracture screw diameter
- Can also be used for other company screw fractures
- Recommended tightening torque : 35Ncm
- Recommended speed : 100rpm or less



| D (Removal Screw) | Ø 1.2 | Ø 1.4 | Ø 1.6 | Ø 1.8 | Ø 2.0 |
|-------------------|----------|----------|----------|----------|----------|
| | OSRT 12E | OSRT 14E | OSRT 16E | OSRT 18E | OSRT 20E |

Wireless Electric Driver

e-Driver 06.2016

- Strong and accurate torque (5~35Ncm)
- Adjustable rpm (15~60rpm)
- Minimizes OrthAnchor fractures and accurate insertion path
- Easy abutment tightening and minimizing the chance of screw loosening

OSM-TORQ



e-Driver Plus NEW 03.2021

- Easy tightening of contra angles
- Strong and accurate torque (5~40Ncm)
- Adjustable rpm (15~55rpm)
- Minimizes OrthAnchor fractures and accurate placement path
- Easy abutment tightening and minimizing the chance of screw loosening


DSD-DTD-0100



V-ceph

V-ceph

- Orthodontics diagnostic software
- VTO / STO (simulation of facial changes before and after treatment)
- Grid view (check the symmetry in the front picture with the guide line)
- Dual monitor-views (compare patient data on two monitors)
- X-ray superimposition (with tracing)
- Sticky note (note in all image views)
- Gallery format (23 types)
- Image process (Image editing)
- Growth forecast
- Change axis (fix FH line horizontally)
- Smart V-ceph (iPad application)

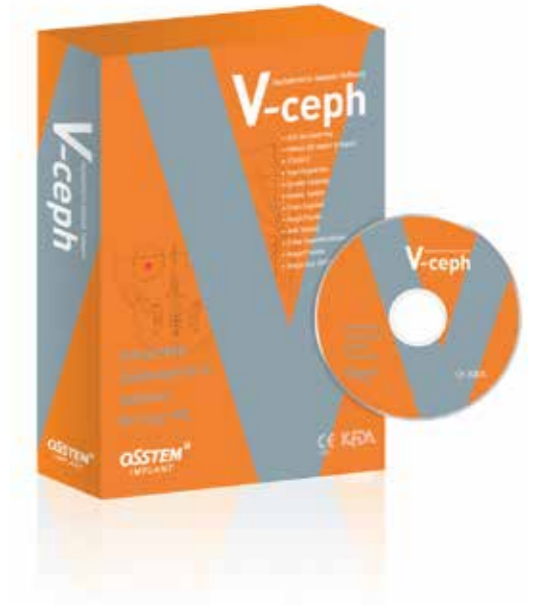

OSSTEM[®]
IMPLANT